These statements have not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure or prevent any disease. This device does not cure cancer. It does not intend to provide diagnosis or treatment and only claims the statements in the medical device listing of intended use approved by the FDA. Specific medical advice should be obtained from a licensed health care practitioner. The information and personal testimonies about this device do not reflect any medical claims.

RICHWAY The name you can trust!

We believe that people must be healed from their wounds, renewed in spirit, restored from illness, enlightened from ignorance, and saved from pain. No one should be left behind. Our beauty does not lie in our appearance, but in our passion of caring for others.

THE RICHWAY GROUP

A State of The Art

AMETHYST BIOMAT 7000MX

THE BIOMAT IS A PRODUCT OF MODERN TECHNOLOGY WITH PROVEN EFFICACY CRAFTED BY HIGHLY SKILLED SCIENTISTS AND ENGINEERS

A BREAKTHROUGH



AMETHYST BIOMAT

BIOMAT 7000mx

Natural Amethyst Biomat US FDA Medical Device Japan FDA Medical Device Korea FDA Medical Device Europe Medical Device CE for Europe / CU for Russia CQC for China


RICHWAY

RICHWAY & FUJI BIO INC.

1314 South King Street Suite 520, Honolulu, HI 96814 USA Tel: 808-589-2800 Toll-Free: 855-338-6410 Fax: 808-597-1651 www.richwayandfujibio.com www.richwaykorea.com



Richway Life with a Biomat

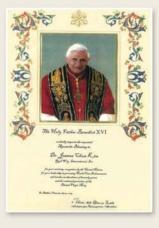

Pope Benedict XVI's blessing for Richway's contribution to health care enhancement.

"The Holy Father Benedict XVI cordially imparts the Apostolic Blessing to Dr. Jeanna Chun Kim, Richway International, Inc., for your receiving recognition by the United Nations, for your leadership in promoting Health Care Enhancement and invokes an abundance of heavenly graces and the continued protection of the Blessed Virgin Mary."

May 6, 2009, The Vatican

THE UNITED NATIONS ASSOCIATION GITATION OF SPECIAL MERIT Dr Jeanna Chun Kim, C.E.O. RichWay International, Inc.

RICHWAY

In witness thereof in recognition whereof, I have keresuto put my hand and affixed the Official Seul of the United Nations Association

In	
Put forest Research THE PLAN	
photo, but series	
EMINED NETRONY ANSOLD	ATOM BICADQUARTERS
PRESS Policie des Nations Bierre (24.) s El 11 General (10, Socialization)	HELDA Luine Valuer, Place Boom Di 1-1177 New Sock, NY 2001, UNA

The United Nations Association Citation of Special Merit

Dr. Jeanna Chun Kim, C.E.O. Richway International, Inc.

"In recognition of her outstanding contributions, dedication in promoting people's health. She has also demonstrated strong leadership and commitments to Primary Health Care Enhancement for the benefit of the needy people and the society at large. In witness thereof in recognition whereof, I have hereunto put my hand and affixed the Official Seal of the United Nations Association, on the 16th of June, 2009, in Jakarta, Indonesia."

Prof. Kernal Roemawi MD, Ph.D.



Recognition from Prince Fushimi of the Japanese Cultural Agency

Richway Japan received an award from Prince Fushimi of the Japanese Cultural Agency for its contribution to promoting health care for the Japanese people.



Richway CEO Jeanna Kim and the Queen of Malaysia The Malaysian Royal Family uses the Biomat.



Richway's invitation by the U.S. White House In 2006, Richway was recognized as an excellent company by the President of the United States and the founder was invited to the White House for a dinner party hosted by President Bush.

Richway

Spreading the Message around the World

GREETINGS FROM RICHWAY IN HAWAII

Aloha! Thank you for your recent interest in Richway & Fuji Bio Inc.

We hope that you will find our products of exceptional quality and receive many health benefits from its use. Richway & Fuji Bio Inc. is always seeking to enhance and innovate ways that allow you to live and look both healthy and happy. We create all of our products with cutting-edge technology, engineering, and methodology to meet your expectations of quality, performance and professional design. The Biomat is used by thousands of people all over the world. It has helped to make satisfied users healthier and we feel confident that it will help you with your health concerns as well. Again, thank you for your interest. We appreciate your patronage and wish you a long healthy life with many "rich" rewards.

BRAZII

Mahalo! Many Warm Regards,

EOUATO

Richway & Fuji Bio Inc., Customer Service



Dear Richway Family

We believe in the power of teamwork. We believe that with hard work and passion, anything is possible. We strive for excellence and integrity in everything we do. We recognize the relationship between our lifestyle and our health, and the importance of maintaining wellness.

You want to better yourself. You have a burning desire to succeed. You demonstrate the most important factors for success: commitment, a passion to make a difference, a vision for achieving positive change, and the courage to take action.

We know our purpose is to make life better and we hold to our basic creed: look better, feel better, and live better. Our products reflect the scientific advancements of today. We work to achieve freedom of time and finances. We serve our communities, we respect our environment, and we place our families first.

We strive to be loving, caring human beings, and we teach what we've learned to anyone who is open to it. We believe in Richway, we believe in what we offer, and we believe in our purpose.

Warmest regards,

Chen C. Am

Jeanna Kim President & CEO Richway & Fuji Bio Inc.



biomat 7000mx

UNIQUE INNOVATION IN BIO TECHNOLOGY

The Biomat's 17 layers consist of:

- 1 Surface material: Silicon urethane with cotton
- 2 Waterproof layer
- 3 Amethyst layer for transferring natural infrared rays
- Hyron cotton layer for thermal insulation
- **5** TOCA layer for natural negative ions
- Nano Copper fabric layer
- Quantum Energy layer (peach and grape seeds)
- (B) Copper fabric layer for electromagnetic interception
- (9) Carbon fiber layer for electromagnetic interception
- Fiberglass layer
- 1 Thermal preservation layer
- **1** Silicon and Teflon reverse currency heating layer with EMF interception
- 13 Nonwoven fabric layer
- (I) Aluminum layer for reflection of infrared rays
- ID Nonwoven fabric layer for heat preservation
- (1) Thermal protection layer
- 1 Bottom material: High quality cotton with brass pattern



SAFETY SYSTEM



CONTROLLER FEATURES

SOFT TOUCH POWER SWITCH AND TEMPERATURE/TIMER CONTROL

The Biomat 7000mx's control system has many useful functions for setting the temperature, timer, and it even has the most advanced I.C. chip to reduce electromagnetic waves.

F.I.R. Indicator (Far Infrared Ray)

This light indicates when far infrared rays are electrically produced. This light will come on when you first set the Biomat's temperature. Once the Biomat has reached the set temperature, this light will turn off and the • E.D.N.I. light will begin to flash.

E.M.F. Interception Light

This lights up to indicate the active reduction of harmful electromagnetic waves. It turns off if there are no electromagnetic waves to block.

Alarm Light

This light will appear if there is any malfunction. If so, check connections to ensure they are secure. Contact the dealer if you are experiencing any issues.

Temperature

Indicator

(Celsius)

Temperature Up

Temperature • Down

Soft Touch
Power Switch



Awarded Best Medical Device in 2010

Grand prize in the medical device field by the Royal Swedish Academy of Sciences





Japan FDA Medical Device BG 10300403

Korea MFDS Medical Device A83080.01[2]

Temperature Indicator (Fahrenheit)

Timer Control

This indicates the usage time before the Biomat turns off. Pressing the button will change the timer setting.





The Nobel Prize in Physiology or Medicine 1991

FOR THEIR DISCOVERIES CONCERNING THE FUNCTION OF SINGLE ION CHANNELS IN CELLS



fur Biophysikalische Chemie Gottingen,



Max-Planck-Institut fur Medizinische Fors-chung Heidelberg,

Bert Sakmann, Germany

Press Release: The 1991 Nobel Prize in Physiology or Medicine NOBELFORSAMLINGEN KAROLINSKA INSTITUTE THE NOBEL ASSEMBLY AT THE KAROLINSKA INSTITUTE



The Nobel Assembly at the Karolinska Institute has decided to award the Nobel Prize in Physiology or Medicine for 1991 jointly to Erwin Neher and Bert Sakmann for their discoveries concerning "the function of single ion channels in cells". According to The Cell, their research established the existence of ion channels in cell membranes, some of which permit the flow of positively and negatively charged ions, leading to the examination of a broad range of cellular functions. Eventually this led to discovering the role that ion channels play in diseases and certain disorders, enabling the development of new therapies.

For more information: http://www.nobelprize.org/nobel_prizes/medicine/laureates/1991/

Richway Biomat Certificates






TÜV's EMC Safety Approval

US Trademark Registration Number 2555064

Product testing is about testing tomorrow's products today to ensure they comply with regulatory requirements and industry standards. It also identifies how the manufacturer's product will perform in the marketplace for the end user or consumer. TÜV also compares similar products to ensure the manufacturer's product has a competitive edge. Testing helps companies to develop safer, greener products and to run safer operations with better integrity and a lower social and environmental impact. TÜV helps companies to comply with the safety, quality and environmental regulations of today and in the future. TÜV helps companies build trust and value in their businesses!


Comparison of Normal Heat and Biomat Far Infrared Rays

After using a conventional heating pad • for 10 minutes Only the outer surface of the skin is •

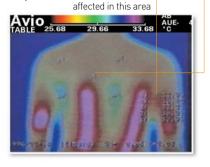

• After 5 minutes of using far infrared rays

Far infrared rays can penetrate the skin much deeper compared to other types of heat

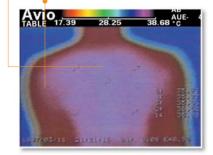

There are two types of heating applications here. One is the simple application of a conventional heating pad, while the other is the result of far infrared rays. One of the main focuses of attention during the 21st century will be on these methods of heating the body using far infrared rays.

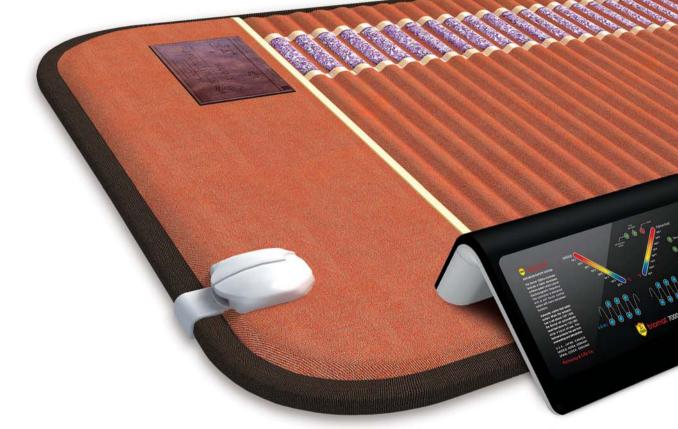


Image: Second Second



U.S. FOA U.S. Food and Drug Administration Approved Medical Device



KOREA MFDS Korea Ministry of Food and Drug Safety Medical Device A83080.01[2]



JAPAN FDA Japan Ministry of Health, Labour and Welfare Medical Device BG10300403

RAL Co. Ltd.
14 An American States Test & M. Sarra Antoproves 1. Same and the approximation of the states of the second states of the states of the states of the states.
performance of the state of the
EN ISO 13485:2012
Note that the second se
PLANE THE REAL PROPERTY

ISO 13485 International Medical Device Manufacturer



ISO 9001 International Electric Equipment Manufacturer

L	100000- 475-440 10000	7-62307
l	의료가가 제도 및 총달관리 ICertificate of 1	718 481984 Call
l	• 200 100 10000 10000 10000 10000 000, 30, 100 • 201 100000000	- @
l	00 000 000 • 10 000 000 000 00000 E10 000 000 000 00000 000 0000 0000 0000	
l	·	
I		
l	· · · · · · · · · · · · · · · · · · ·	-
I	And all a ground a stand of	Not a set of the set of the set
l	KEDA	
l	KTR BARREN	1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1


GMP Good Manufacturing Practice Approval Number 3965 (KGMP : 1-D2337)

What is EMF?

EMF (Electromagnetic Fields) is a term used to describe the invisible electric and magnetic fields radiating away from manmade alternating electrical current. EMF radiates from all electrical appliances in the home, office or factory, in cars, trucks, buses and overhead tram and train lines.

What does the Biomat EMF interceptor do?


The Biomat EMF Interceptor acts in the same way as a lightning rod by safely conducting the EMF energy away to the ground. With this feature, you can effectively reduce exposure to EMF for yourself and your family.

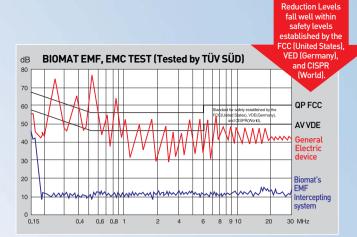


Biomat EMF

These images show circulation before and after 14 days of Biomat usage. Note the increase in circulation.

According to the American Heart Association, the circulatory system is responsible for the delivery of oxygen-rich blood and the removal of wastes. In Dr. Inoue and Dr. Kabaya's paper entitled "Biological activities caused by far-infrared radiation", far infrared rays have been found to have a blood circulation-enhancing effect in human skin and eventually induces an increase in temperature of the body tissues.

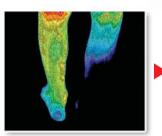

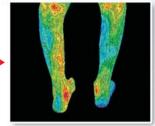



International Electrical Safety Approval Number EUI-SI CE2013-0142



Korean Electrical Safety Korean Electrical Safety Number PRO: HH071485-13002A MINI: HH071485-13001A

IEC IRE	la fan an a
Part 2 17 Parts	Tabl Appoint at MCAL-11 c which exciting application ibolicy der requirements for blankack, posits and in function treating application
Paget Rates	NAME OF TAXABLE PARTY.
Particular in American Internet	La has interested
Papersonal Ra 1 Augustures	it in Particia
Pade of Sector	ana.a
The other of page	Thi cape
Patheolis (and	No.010
-	the party of the party of the party of the
Two and designs	
Second	the second of the second secon
	THE ADDRESS SECTOR
The pinister	
	-
Subardon b	#1808.1.10
The Agent Surger Without	-
man ter	Sectors.
AND DESCRIPTION OF	 A second being and inclusion of the matching in the second of the second second second se second second sec
	and the residues the BORDER's against the effective to the D
The sport over shift are by the	the other states and the second states and t
"will belt impringibut	
Tel bit meriphe	Page land
	en su les

Japan Electrical Safety Japan Electrical Safety Number HW 2013-0145

🕸 ВЕТТ	ER BUSINESS BUREAU
BBB serving Henryl	None Concerner (First Concerner Con
	erving Hawaii Reliability Report
	Extremy International, Inc.
General Information	KOTWAY INSTANDARI, INC.
Need Address.	1314 South King Street Sole 520 Resolute RD 99924 U.S.C.
Kdiphone Number:	00079-200
Kaul Address	Advey was not
Web Address	Next we we have a parallely black as an
Product Officerals	This sempery effect almostic facial & holp care health & beauty servi-
Citra Contractor	
Lines	General Decise License #092197#
Data Original Davisore Statust	Ady pixt
Resid Calmer Control	Mr. Cahar Star Parallel
Program Participation	als Microbenialy Manufacture Program, association franching procedures, including mediation and arbitraries, if necessary, so
	is stillable y solid with the Research Regime Improved encodence on the Research company much have been information for at last (2 modify, pro- by the Research of the One Research measured reference experision of comparison physics excluded, to addition, the Research measurements of the origination lastic.
Licensing Information The sense of the sense of the sense The sense sense is a sense of the sens	iniar Linna An An Annahm at an Inn an An Annahm All Spennarson (von ofen son

BBB Credit Rating(RICHWAY2011)

EAC	SEPTIMOLIAT CON	
		State Bring St.
Income of the local division of the local di		Contraction of the second
for delay	And I want to be a surrow of	Construction of the second second
Statement and	and the second second second	
	A STATISTICS	
1000		
		The second second
And the Party of t		
and .		
		all at here and

Russia CU Electrical Safety License



Antibacterial Activity Safety Test Industry & Environment Research Center

RICHWAY biomot 7000 A STATE OF THE ART

Amethyst Biomat Professional

Includes: Biomat, Controller, Cotton Pad, Carry Case



Amethyst BMP Size Net Weight Electric Consumption Voltage Performance Amethyst Size Amethyst Size Amethyst Weight Amethyst Color Cut & Shape Polariscope Test Refractive Index 700 x 1850 mm / 27.56" x 72.83" 11.11 kg / 24.49 lbs. without controller 180 W AC 120 Temperature 35-70°C / 95-158°F 5 mm - 12 mm 8 kg / 17.63 lbs. Natural, Violet, Clear Tumbled D.R. 1.544-1.553

Amethyst Single Biomat

Includes: Biomat, Controller, Cotton Pad, Bag



Amethyst Single Size Net Weight Electric Consumption Voltage Performance Amethyst Size Amethyst Size Amethyst Weight Amethyst Color Cut & Shape Polariscope Test Refractive Index

Amethyst Queen Size

Net Weight

1000 x 2000 mm / 39.37" x 78.74" 17.18 kg / 37.87 lbs. without controller 180 W AC 120 Temperature 35-70°C / 95-158°F 5 mm - 12 mm 9 kg / 19.84 lbs. Natural, Violet, Clear Tumbled D.R. 1.544-1.553

Amethyst Queen Biomat

Includes: Biomat, Controller, Cotton Pad, Bag



Electric Consumption Voltage Performance Amethyst Size Amethyst Weight Amethyst Color Cut & Shape Polariscope Test Refractive Index

Amethyst King Size Net Weight Electric Consumption Voltage Performance Amethyst Size Amethyst Size Amethyst Weight Amethyst Color Cut & Shape Polariscope Test Refractive Index 1400 x 2000 mm / 55.12" x 78.74" 27.1 kg / 59.74 lbs. without controller 360 W AC 120 Temperature 35-70°C / 95-158°F 5 mm - 12 mm 14 kg / 30.86 lbs. Natural, Violet, Clear Tumbled D.R. 1.544-1.553

1800 x 2000 mm / 70.87" x 78.74" 33.33 kg / 73.48 lbs. without controller 360 W AC 120 Temperature 35-70°C / 95-158°F 5 mm - 12 mm 20 kg / 44.09 lbs. Natural, Violet, Clear Tumbled D.R. 1.544-1.553

Amethyst King Biomat

Includes: Biomat, Controller, Cotton Pad, Bag



A BREAKTHROUGH IN MEDICAL SCIENCE

Amethyst & Tourmaline Mini Biomat

Includes: Biomat, Controller, Cotton Pad, Carry Case



Amethyst & Tourmaline Pillow

Includes: Memory Foam Pillow, Amethyst & Tourmaline Cover, Cotton Cover





Amethyst Mini Size Net Weight Electric Consumption Voltage Performance Amethyst Size Amethyst Veight Amethyst Color Tourmaline Size Tourmaline Weight Tourmaline Color Cut & Shape Polariscope Test Refractive Index 500 x 850 mm / 19.69" x 33.46" 3.6 kg / 7.93 lbs. without controller 100 W AC 120 Temperature 40-70°C / 104-158°F 5 mm - 12 mm 1.2 kg / 2.64 lbs. Natural, Violet, Clear 5 mm - 12 mm 1.2 kg / 2.64 lbs. Black Tumbled D.R. 1.544-1.553

Pillow Size Net Weight Amethyst Size Amethyst Weight Amethyst Color Tourmaline Size Tourmaline Weight Tourmaline Color 480 x 300 x 120 mm / 19" x 12" x 4.3" 2.5 kg / 5.51 lbs. 5 mm - 12 mm 0.6 kg / 1.32 lbs. Natural, Violet, Clear 5 mm - 12 mm 0.6 kg / 1.32 lbs. Black

Elastic Polyurethane Memory Foam Pillow

Tourmaline has a complex chemical composition, abundant in boron and aluminum silicate minerals of variable composition. Tourmaline comes in a wide range of colors, depending on its mineral composition.

Size	Comforter Pillow cover	2000 x 2300 mm / 78.74" x 90.55" 790 x 500 mm / 31.10" x 19.69"
Тор	100% Polyester Jacquard	
Bottom	100% Spining	
Inner	100% Polyeste	er (Microfiber)

Microfiber: The ultimate product of nanofiber technology!

Richway's development of the Quantum Energy Comforter, with its combination of warm and light microfiber and advanced quantum energy technology, will give you the most comfortable night's sleep, especially combined with the Biomat.

Size	Pro 700 x 1850 mm / 27.56" x 72.83"	
	Single 1000 x 2000 mm / 39.37" x 78.74"	
	Queen 1500 x 2000 mm / 59.06" x 78.74"	
	King 1800 x 2000 mm / 70.87" x 78.74"	
Cover	Microfiber Fabric, 100% Polyester	
Cushion	100% Polyester, 3D Mesh	

3D Mesh is made of a dense, regular molecular structure with outstanding air permeability. The comfortable 3D Mesh distributes body weight evenly and removes pressure points, providing luxurious comfort. It also works great for those who own a Biomat and prefer more cushion support.